

Production: LEMON CURD

Objective: Comparison of spreadability/firmness of luxury and economy lemon curd

Type of action: Penetration test

Test setting:

| Speed | Test mode | Trigger | Target | Hold |
|--------|--------------|---------|--------|-------|
| 3 mm/s | Distance (c) | 0 gf | 23 mm | 0 sec |

Accessory:

Butter spreadability rig, Platform

Accessory:

Many people are interested in measuring the spreadability of food products such as butter, spreads, peanut butter, margarines, cheese, and cream cheese as well as non-food items such as lotions, creams, gel shaving creams, joint compound, greases, etc. These products are often very shear sensitive and are difficult to consistently prepare for testing. Penetration and compression style tests are simple methods that give results of sample hardness. Hardness measurements however, even with cone probes, are not always good discriminators of spreadability although the apparent yield value calculated from cone penetration with constant load has been used in the past and has correlated with sensory estimates of spreadability.

Firmness (or hardness - the two are synonymous) as a discrete sensory property of a tablespread is implicit in the determination of its spreadability. The first impression of a tablespread texture will be that of its firmness, obtained as the consumer forces a spreading instrument into the product. Spreadability, in pragmatic terms, is the ease with which a spread can be applied in a thin, even layer to bread. Firmness or hardness may be measured by the force required to obtain a given deformation or by the amount of deformation under a given force. Although spreadability is also a deformation under an external load, it is a more dynamic property. Measurements of firmness and spreadability are usually highly correlated, however the relationship is rarely perfect, and this is partly a function of work-softening. Margarine, for instance, worksoftens (when spread on bread) more easily than butter, which allows it to be more spreadable even when hardness values are initially equal.

The Butter spreadability rig is a set of precisely matched male and female perspex 90 cones. The material is allowed to set up in the lower cone or is filled into the lower cone with a spatula. The material is pressed down only so much as is needed to eliminate air pockets which are visible through the perspex cones, and then the surface is levelled with a flat knife. Excessive work is not introduced into the product. The fixture comes with five replaceable female cone sample holders which can be filled in advance of testing and then easily locked into the base holder precisely centred under the matching upper male cone probe. The sample holders can be stored in refrigerated or frozen environments, or they can be used at room temperature.

As long as the material has not been excessively worked or whipped, different styles of filling the material only affect the early part of the test.

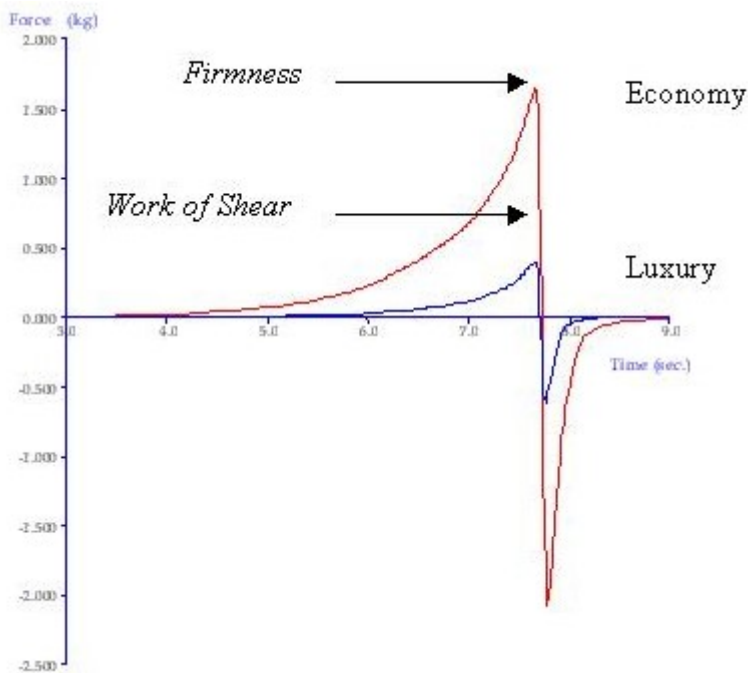
Sample Preparation:

Place sample into the female cone and press it down to eliminate air pockets. Before filling all of the sample holders, use an empty holder to perform the probe calibration procedure (see below). Scrape off any excess sample with a knife, to leave a flat test area. Excessive work should not be introduced into the product. Allow the samples to equilibrate to the specified temperature e.g. 5C before testing.

Test Set-Up:

Place the Heavy Duty Platform on the base of the machine but do not secure the position at this stage. Position the base holder on the Heavy Duty Platform and lock into place by tightening the screws. Insert an empty female cone sample holder into the base holder. Move the male cone probe down so that it fits almost all of the way into the female cone sample holder. When the male and female cones are practically touching, manoeuvre the Heavy Duty Platform so that the cones are precisely aligned. At this point tighten the screws of the Heavy Duty Platform so that it is secured into position. Before testing, the male cone probe must be calibrated against the female cone so that the starting point is at the same height for each test e.g 25.0mm above the female cone. To do this, raise the male cone up so that it is just above the female cone and select: **Calibrate Height**.

When prompted, type in 25.0mm for the return distance, then click on OK. The probe will come down, just touch the female cone, and then return to a position precisely 25.0mm over the female cone. The **Control Probe** facility allows this test position to be reset rapidly and accurately after replacing the female cone with a new sample and cleaning

Typical plots:

The above curves were produced from luxury and economy lemon curd tested at 5C.

Observations:

Once the test is commenced the male cone probe proceeds to approach and then penetrate into the sample and continues to a depth of 2mm above the sample holder surfaces, i.e. the probe moves a distance of 23mm from its start point. During penetration the force is seen to increase up until the point of maximum penetration depth. This force value can be taken as the "Firmness" at this specified depth. A firmer sample also shows a correspondingly larger area that represents the total amount of force required to perform the shearing process. Both of these values have been shown to rank samples in the same order of spreadability (and firmness), but for some samples one may prove to be more suitable than the other. The probe then proceeds to withdraw from the sample.

Data Analysis:☒ Max Force☒ Area (+)**Results**

| Sample | Mean Max. +ve Force 'Firmness' (+/- S.D.) (kg) | Mean Area 'Work of Shear' (+/- S.D.) (kg s) |
|---------|---|--|
| Economy | 1.6 +/- 0.2 | 1.4 +/- 0.2 |
| Luxury | 0.4 +/- 0.03 | 0.2 +/- 0.01 |

Notes:

- To avoid the Texture Analyser from overloading/underloading, the cones must be precisely aligned and the chosen penetration distance should always be at least two millimetres less than the calibrated probe height.
- In between tests, only remove the female cone from the base holder; do not remove the base holder from the Heavy Duty Platform otherwise the male and female cones will require realignment.
- During penetration, a large blip in the otherwise smooth curve may be observed. This is due to the compression of an air pocket within the product, so it is important to minimise the presence of air pockets when filling the female cone.
- When reporting results the test temperature and calibrated height (start point) should always be specified for the purpose of comparison.
- When attempting to optimise test settings it is suggested that the first tests are performed on the hardest samples to anticipate the maximum testing range required and ensure that the force capacity allows testing of all future samples.